




USER MANUAL – EN

**IN 27604 Facial LED light tunnel inSPORTline Coladome
300**

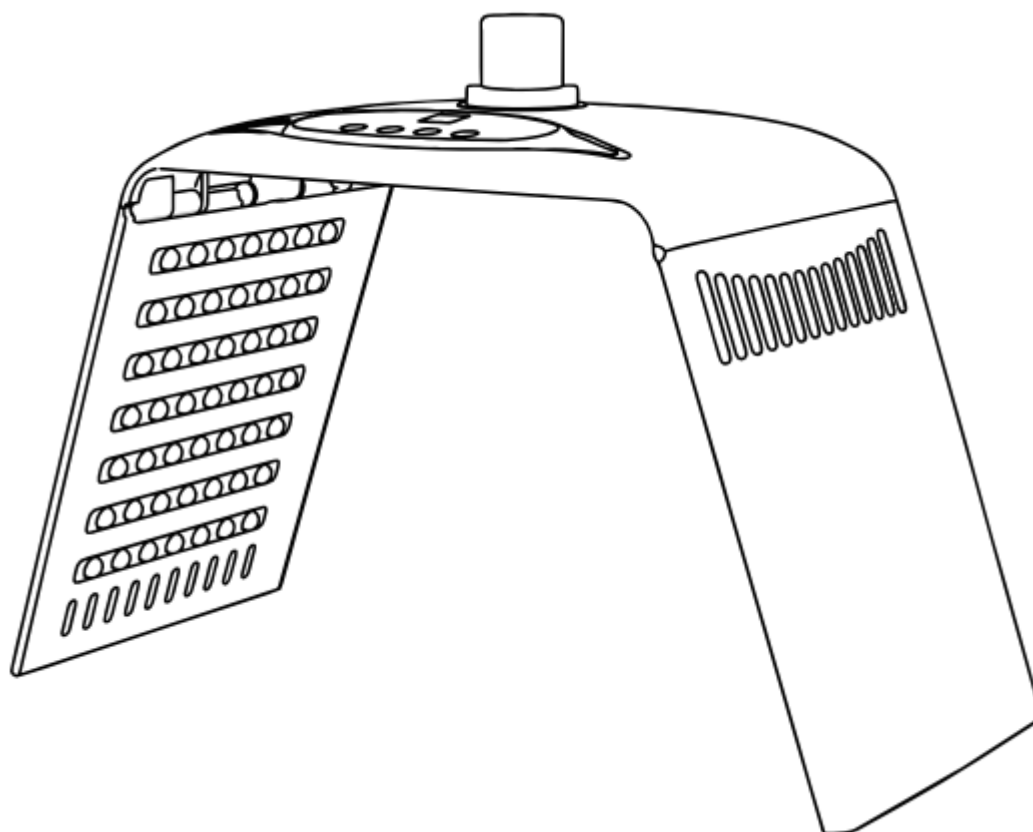

CONTENTS

SAFETY INSTRUCTIONS..... 3
PRODUCT DESCRIPTION 4
UNFOLDING..... 4
CONTROL PANEL 5
SPECIFICATIONS..... 6
MAINTENANCE 6
ENVIRONMENT PROTECTION 6
TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS..... 6

SevenSport s.r.o. reserves the right to make any changes and improvements to its product without prior notice. Visit our website www.insportline.eu where you will find the latest version of the manual.

SAFETY INSTRUCTIONS

- Read the manual before use and keep it for future reference.
- Keep the manual for other users to view.
- Inform other users about basic safety instructions and use.
- The device may only be used by an experienced and adult person with required knowledge.
- Prior to use, make sure that the LED lights and the device are not damaged.
- Do not use or store the device in places with high humidity (bathrooms, saunas, swimming pools, etc.).
- Use only in a well-ventilated room.
- The device must be placed on a flat, firm, and stable surface.
- Protect the device from impacts or falls.
- Do not use other heat radiating devices during use.
- Keep out of reach of water.
- Do not plug into a damaged outlet.
- Do not use if the cord or plug are damaged.
- After use, switch off the device and unplug it. Unplug the device for maintenance and cleaning.
- If water enters the device, unplug it immediately from current.
- If the device makes abnormal noises or smells, stop operation immediately and unplug the device from current.
- Do not modify or repair the device yourself.
- Use the device only as described in the manual. Any other use is not permitted.
- Do not place the device near heat sources or open fire.
- Do not cover with a cloth or other objects during use.
- Only clean the device with a dry cloth.
- Do not use if you suffer from skin diseases, heart problems, high blood pressure, irritation, or sensitive skin.
- Always wear safety goggles!
- Recommended operating temperature: 0 - 26 ° C
- Recommended storage temperature 0 - 35 ° C
- Relative humidity: <80%
- Frequency of use: can be used 25-30 minutes a day for one month or after any cosmetic procedure. As part of acne treatment, we recommend using it 3 times a week, for at least 25 minutes

PRODUCT DESCRIPTION



1. Control panel	5. LEDs
2. Water container	6. Power cable
3. Power supply	7. Safety glasses
4. Steam outlet	


UNFOLDING

Before unfolding, it is necessary to select a suitable place away from heat sources and other devices. Keep away from sunlight and high humidity. Use in a well-ventilated room with stable temperature and humidity.



1. Unfold both sides of the machine and place on a flat and firm surface.
2. Plug the power cord into the device.
3. Plug the power cord into the electrical outlet.
4. After switching on, the display will light up.

CONTROL PANEL

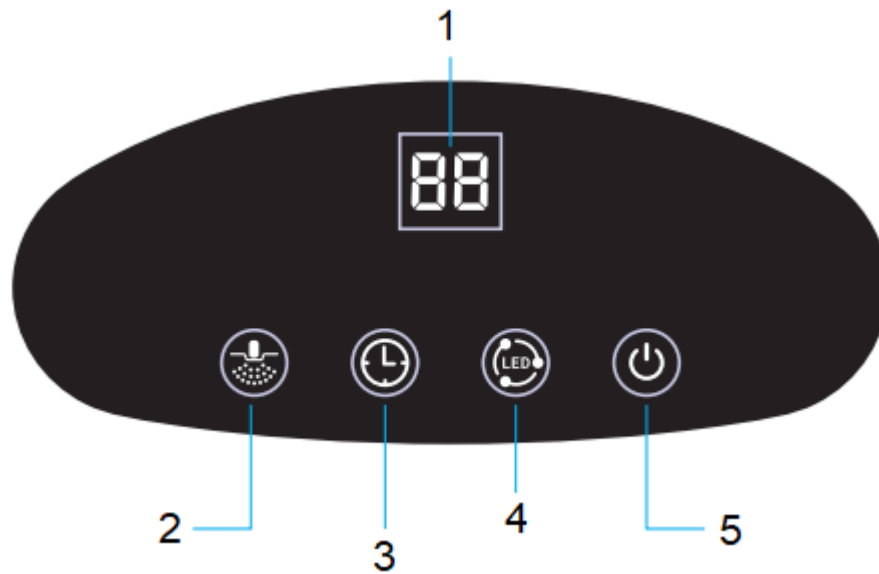

1. Display	4. LED lights
2. Steam function	5. Power switch
3. Timer	

SWITCHING ON

After connecting to the power supply, the device is in standby mode. You start the device by pressing the power switch. The LED indicator and the 30 min timer will light up.

To turn off, press the power button again.

LED

By pressing the button to set the LED lights, you can select one of 8 options.

00 – standby mode

01 – red light (wavelength 647 – 760 nm)

02 – blue light (wavelength 430 – 455 nm)

03 – green light (wavelength 492 – 550 nm)

04 – yellow light (wavelength 550 – 580 nm)

05 – orange light (wavelength 588 – 647 nm)

06 – violet light (wavelength 360 – 430 nm)

07 – teal light (wavelength 480 – 490 nm)

08 – cycling light

TIMER

The initial value of the timer is 30 min. By pressing the button, you can adjust the value by 5 minutes up to a maximum value of 60 min.

Time setting range 0 – 60 minutes.

If the value is not adjusted for 3 seconds, the device will start working on the currently set value.



STEAM FUNCTION

After switching on, the steam function is turned off. Press the button to switch between intensity (01 lowest, 03 highest).

Do not pour essential oils or other cosmetic products into the container. We recommend using clean water.

SPECIFICATIONS

Power supply	DC12V, 3A, 50/60 Hz
Power	Ca. 36W
Dimensions	470x220x345 mm
Weight	1,3 ±10%

MAINTENANCE

- Wipe the body of the device with a soft cloth.
- Store in a dry and shady place.
- Fold the device for storage.
- Before any handling or storage, the device must be disconnected from the electricity supply.

ENVIRONMENT PROTECTION

After the product lifespan expired or if the possible repairing is uneconomic, dispose it according to the local laws and environmentally friendly in the nearest scrapyards.

By proper disposal you will protect the environment and natural sources. Moreover, you can help protect human health. If you are not sure in correct disposing, ask local authorities to avoid law violation or sanctions.

Don't put the batteries among house waste but hand them in to the recycling place.

TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS

General Conditions of Warranty and Definition of Terms

All Warranty Conditions stated here under determine Warranty Coverage and Warranty Claim Procedure. Conditions of Warranty and Warranty Claims are governed by Act No. 89/2012 Coll. Civil Code, and Act No. 634/1992 Coll., Consumer Protection, as amended, also in cases that are not specified by these Warranty rules.

The seller is SEVEN SPORT s.r.o. with its registered office in Strakonická street 1151/2c, Prague 150 00, Company Registration Number: 26847264, registered in the Trade Register at Regional Court in Prague, Section C, Insert No. 116888.

According to valid legal regulations it depends whether the Buyer is the End Customer or not.

"The Buyer who is the End Customer" or simply the "End Customer" is the legal entity that does not conclude and execute the Contract in order to run or promote his own trade or business activities.

"The Buyer who is not the End Customer" is a Businessman that buys Goods or uses services for the purpose of using the Goods or services for his own business activities. The Buyer conforms to the General Purchase Agreement and business conditions.

These Conditions of Warranty and Warranty Claims are an integral part of every Purchase Agreement made between the Seller and the Buyer. All Warranty Conditions are valid and binding, unless

otherwise specified in the Purchase Agreement, in the Amendment to this Contract or in another written agreement.

Warranty Conditions

Warranty Period

The Seller provides the Buyer a 24 months Warranty for Goods Quality, unless otherwise specified in the Certificate of Warranty, Invoice, Bill of Delivery or other documents related to the Goods. The legal warranty period provided to the Consumer is not affected.

By the Warranty for Goods Quality, the Seller guarantees that the delivered Goods shall be, for a certain period of time, suitable for regular or contracted use, and that the Goods shall maintain its regular or contracted features.

The Warranty does not cover defects resulting from (if applicable):

- User's fault, i.e. product damage caused by unqualified repair work, improper assembly, insufficient insertion of seat post into frame, insufficient tightening of pedals and cranks
- Improper maintenance
- Mechanical damages
- Regular use (e.g. wearing out of rubber and plastic parts, moving mechanisms, joints, wear of brake pads/blocks, chain, tires, cassette/multi wheel etc.)
- Unavoidable event, natural disaster
- Adjustments made by unqualified person
- Improper maintenance, improper placement, damages caused by low or high temperature, water, inappropriate pressure, shocks, intentional changes in design or construction etc.

Warranty Claim Procedure

The Buyer is obliged to check the Goods delivered by the Seller immediately after taking the responsibility for the Goods and its damages, i.e. immediately after its delivery. The Buyer must check the Goods so that he discovers all the defects that can be discovered by such check.

When making a Warranty Claim the Buyer is obliged, on request of the Seller, to prove the purchase and validity of the claim by the Invoice or Bill of Delivery that includes the product's serial number, or eventually by the documents without the serial number. If the Buyer does not prove the validity of the Warranty Claim by these documents, the Seller has the right to reject the Warranty Claim.

If the Buyer gives notice of a defect that is not covered by the Warranty (e.g. in the case that the Warranty Conditions were not fulfilled or in the case of reporting the defect by mistake etc.), the Seller is eligible to require a compensation for all the costs arising from the repair. The cost shall be calculated according to the valid price list of services and transport costs.

If the Seller finds out (by testing) that the product is not damaged, the Warranty Claim is not accepted. The Seller reserves the right to claim a compensation for costs arising from the false Warranty Claim.

In case the Buyer makes a claim about the Goods that is legally covered by the Warranty provided by the Seller, the Seller shall fix the reported defects by means of repair or by the exchange of the damaged part or product for a new one. Based on the agreement of the Buyer, the Seller has the right to exchange the defected Goods for a fully compatible Goods of the same or better technical characteristics. The Seller is entitled to choose the form of the Warranty Claim Procedures described in this paragraph.

The Seller shall settle the Warranty Claim within 30 days after the delivery of the defective Goods, unless a longer period has been agreed upon. The day when the repaired or exchanged Goods is handed over to the Buyer is considered to be the day of the Warranty Claim settlement. When the Seller is not able to settle the Warranty Claim within the agreed period due to the specific nature of the Goods defect, he and the Buyer shall make an agreement about an alternative solution. In case such agreement is not made, the Seller is obliged to provide the Buyer with a financial compensation in the form of a refund.

CZ
SEVEN SPORT s.r.o.

Registered Office: Strakonická 1151/2c, Praha 5, 150 00, ČR
Headquarters: Dělnická 957, Vítkov, 749 01
Warranty & Service: Čermenská 486, Vítkov 749 01

CRN: 26847264
VAT ID: CZ26847264
Phone: +420 556 300 970
E-mail: eshop@insportline.cz
reklamace@insportline.cz
servis@insportline.cz
Web: www.inSPORTline.cz

SK
inSPORTline s.r.o.

Headquarters, warranty & service center: Električná 6471,
Trenčín 911 01, SK

CRN: 36311723
VAT ID: SK2020177082
Phone: +421(0)326 526 701
E-mail: objednavky@insportline.sk
reklamacie@insportline.sk
servis@insportline.sk
Web: www.inSPORTline.sk

About shipping

